

Technical Sheet

Sensor Mobile Support

NEW
MÉDI



Ref. : EIM473


Technical data :

Weight: 11 kg
Height: from 0 to 200 cm
footprint:
H.120 cm x 60 x 60 cm

Compatibility :


DR sensors
CR cassettes
Standard cassettes (film)

fall protection

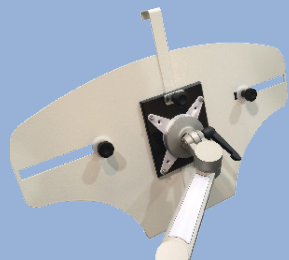

options:


Support for 36x43cm grid
protection. AIM491

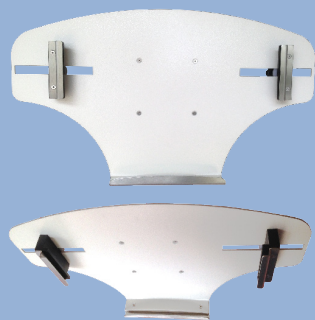

Basket sensor storage (or tape)
to attach to the foot.
AIM492



Capt'X

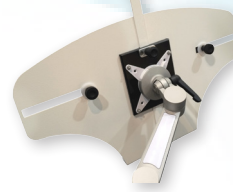

New mobile support «CAPT'X EASY» for DR sensors or CR cassettes, as well as for the standard cassette (film) will get you the position desired of them quickly and easily.

The «CAPT'X EASY» is very versatile and can be used in all services where you will need to do X-rays: emergency, operating room, bedroom, in radio rooms, etc...

The «CAPT'X EASY» is especially suitable for examinations that need to be taken slightly below or above the table top, a stretcher, a wheelchair or bed Very handy for reviews of ankles, feet, knees but also for all other body parts.

Its sub-structure down (foot) gives it great stability, made of steel 6 mm thick with an epoxy finish (stove enamelled). The wheels of a diameter of 100 mm give it great maneuverability.

The vertical movements are telescopic, the system is light enough not to need the jack type of assistance or counterweight, often recognized for mechanical problems generating technician interventions.

His kneecap for positioning the sensor in any desired inclinations. The double arm allowing flexible adjustment, fast and precise depth of the sensor, all this in order to allow each user of the possibilities of optimal configurations.

